


THE OPTIMAL CUSHION FOR WHEELCHAIR USERS WITH BODY ASYMMETRY

Pre-ischial bar limits sliding forward | converts shear into normal forces

Follows the body's contours for equal pressure distribution | support IT-coccyx area | managing pelvic asymmetry

Thigh support

Mid height zip closure for easy adjustment when client is seated | cushion keeps its shape, enhancing the user's skin protection

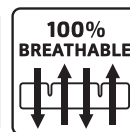

Individual solution for asymmetry

Minimizes sliding and shear forces

Optimal adjustment to the body's contours for equal pressure distribution

The Adjuster O2 is extremely suitable for use with asymmetric body contours and amputations. It automatically manages small obliquities and is individually adjustable to support more extreme asymmetry. When you're looking for a high comfort level, optimal skin protection and stability: the Adjuster O2 is always a safe choice. The design allows total immersion to equally redistribute pressure and limits sliding forward. 100% machine washable and breathable, for optimal hygiene and perfect temperature and moisture regulation. As a result the Adjuster O2 offers you optimal protection against skin breakdown.

PRODUCT PERFORMANCE

|                                  |       |
|----------------------------------|-------|
| Skin protection / Anti decubitus | ▲▲▲▲▲ |
| Pressure redistribution          | ▲▲▲▲▲ |
| Shear force reduction            | ▲▲▲▲▲ |
| Microclimate control             | ▲▲▲▲▲ |
| Impact & vibration damping       | ▲▲▲▲▲ |

Positioning

|                           |        |
|---------------------------|--------|
| Lateral stability         | ▲▲▲▲▲△ |
| Frontal stability         | ▲▲▲▲▲  |
| Sliding reduction         | ▲▲▲▲▲  |
| Allows manual adjustments | ▲▲▲▲▲  |
| Allows easy transfers     | ▲▲▲▲▲△ |

Maintenance & Safety

|                            |        |
|----------------------------|--------|
| Durability                 | ▲▲▲▲▲  |
| Reliability                | ▲▲▲▲▲  |
| Functionality when misused | ▲▲▲▲▲△ |
| Low maintenance level      | ▲▲▲▲▲  |
| Level of hygiene           | ▲▲▲▲▲  |

Functionality & Application

|                          |        |
|--------------------------|--------|
| Comfort                  | ▲▲▲▲▲  |
| Low cushion weight       | ▲▲▲▲▲  |
| Available size range     | ▲▲▲▲▲  |
| Maximum user weight      | ▲▲▲▲▲  |
| Low seat to floor height | ▲▲▲▲▲  |
| Auto adjustments         | ▲▲▲▲▲△ |

USED WITH


Hemiplegia (CVA), amputations, foot propulsion, kyphosis, pelvic obliquity/ tilt/rotation, Spina Bifida, spasms, Fibrodysplasia Ossificans Progressiva (FOP)

TECHNICAL SPECIFICATIONS

|                  |             |
|------------------|-------------|
| Height           | 6 and 10 cm |
| Average weight   | 725 grams   |
| Max. user weight | 250 kg      |
| Machine washable | 60°C max.   |

CUSHION SIZES

Standard (Sizes >53 cm height 10 only)  
On request

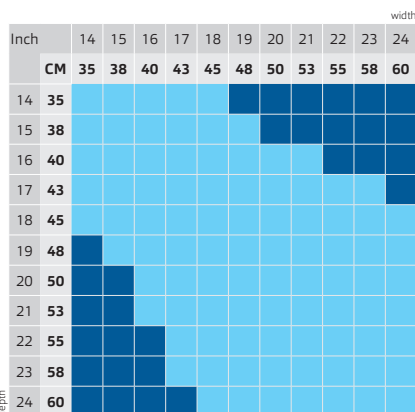

070069-OAD-EN-REV01/2020